
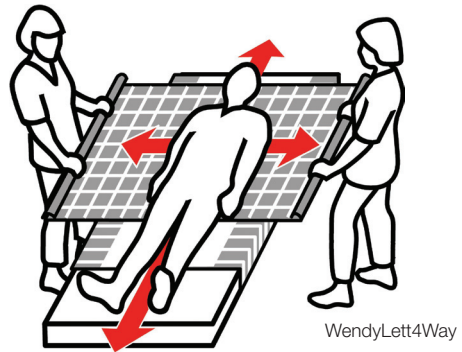


WendyLett2Way

WendyLett4Way

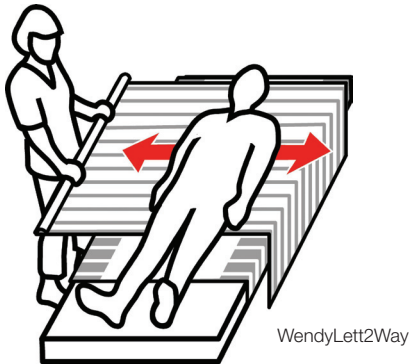

SWL:

300 kg/661 lbs

Art.no. 1649 och 1657:

200 kg/440 lbs

WendyLett2Way and WendyLett4Way are two different draw sheets that can be used together with a WendyLett slide sheet. WendyLett2Way is used for lateral repositioning and WendyLett4Way is used for both lateral repositioning and repositioning higher up or lower down in bed. WendyLett4Way has an anti-friction material on the underside and has been developed mainly for heavy patients. Available in several sizes.

WendyLett is part of a series of different models of transfer assistive devices. The choice of model and material is determined by the patient's functional capacity and the situations in which the product is to be used.



Functional inspection

Visual inspection

Check the condition and function of the assistive device regularly. Always inspect the product after laundering. Check to ensure that seams and material are free from damage. Check to ensure that fabric is not worn or faded. If there are signs of wear, the product must be discarded.




Always read the manual

Always read the manuals for all assistive devices used during a transfer.

Keep the manual where it is accessible to users of the product.


Do not leave the patient unattended during a transfer situation.

Placement

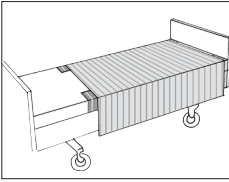

WendyLett2Way

Place WendyLett with the smooth surface in the middle of the bed and with WendyLett2Way on top. After turning, gliding movement can be inhibited by bedding WendyLett2Way in under the mattress on both sides of the bed.

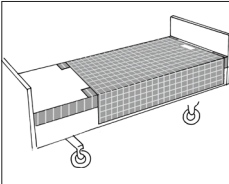

WendyLett4Way

A tag on WendyLett4Way indicates which side should be placed facing upwards in bed. After turning, gliding movement can be inhibited by bedding WendyLett4Way in under the mattress on both sides of the bed.

Turning with WendyLett2Way



Perform the turn by grasping the draw sheet at the level of the hips and the shoulders. For best ergonomic technique, grasp a lot of fabric.



Work with straight wrists and shift weight by leaning backwards.



Turning with the aid of a draw sheet can be performed by two caregivers. This is a collaborative effort whereby one caregiver pulls the draw sheet (as above) and the other presses down somewhat near the shoulders and hips while at the same time pushing.




Turn the patient. After turning, gliding is prevented by tucking the draw sheet under the mattress on both sides.

Turning with WendyLett4Way



The WendyLett4Way draw sheet, which has an anti-friction material on the underside, has been developed specially for heavy patients. The procedure for turning is the same as for WendyLett2Way (above).

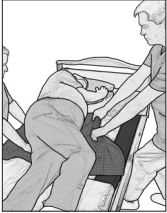

After turning, gliding is prevented by tucking the draw sheet under the mattress on both sides. To prevent chafing, smooth out the bedding under the patient before tucking the draw sheet in under the mattress. The caregivers can do this by pulling diagonally on the draw sheet from either side. Raise the bed rails before leaving the patient.

Higher up in bed with WendyLett4Way



To reposition the patient higher up in bed, grasp the draw sheet on both sides near the shoulders, stand in a walking stance and take one step backwards.



To prevent the patient from gliding down in bed, tuck in part of the draw sheet, so that it is folded double under the patient's seat.

Sitting up in bed with WendyLett4Way



With a WendyLett4Way in the bed, downward gliding can be prevented even when the patient sits up in bed. Tuck in part of the draw sheet, so that it is folded double under the patient's seat. Then, tuck the draw sheet in under the mattress, but only at the hip section. When the head end of the bed is subsequently raised, the draw sheet acts as a friction surface to prevent sliding.



Note:



Be aware that there is a risk of the patient sliding out of bed. Raise the bed rails before leaving the patient.

Material:



The WendyLett slide sheet has a pattern of white stripes on a grey background. Stripes indicate direction of glide.



The WendyLett2Way draw sheet has a white fabric with grey stripes. Stripes indicate the direction in which the sheet is drawn.



WendyLett4Way has a grey-checked pattern and has a sliding surface underneath. Stripes indicate the direction in which the sheet is drawn.

Polyester, cotton, satin

Care of the product:

Read the product label.



Do not use rinsing agent. To ensure maximum material life, hang dry or tumble dry on low heat.

Item No. Product

Size

Colour

WendyLett - sliding sheet

(Can be combined with WendyLett2Way drawsheet or WendyLett4Way sliding drawsheet)

1609	WendyLett, 60 cm/24" sliding surface (<i>bed 80/90 cm - 32"/36"</i>)	100 x 200 cm/40" x 80"	white
1614	WendyLett, 60 cm/24" sliding surface (<i>bed 80/90 cm - 32"/36"</i>)	140 x 200 cm/56" x 80"	white
1629	WendyLett, 60 cm/24" sliding surface (<i>bed 80/90 cm - 32"/36"</i>)	100 x 200 cm/40" x 80"	grey/white striped
1634	WendyLett, 60 cm/24" sliding surface (<i>bed 80/90 cm - 32"/36"</i>)	140 x 200 cm/56" x 80"	grey/white striped
1639	WendyLett, 70 cm/28" sliding surface (<i>bed 105 cm/42"</i>)	200 x 200 cm/80" x 80"	grey/white striped
1667	WendyLett, 100 cm/40" sliding surface (<i>bed 120 cm/48"</i>)	100 x 220 cm/40" x 88"	grey/white striped
1668	WendyLett, 100 cm/40" sliding surface (<i>bed 120 cm/48"</i>)	140 x 220 cm/56" x 88"	grey/white striped
1669	WendyLett, 100 cm/40" sliding surface (<i>bed 120 cm/48"</i>)	200 x 220 cm/80" x 88"	grey/white striped
1641	WendyLett, 70 cm/28" sliding surface (<i>bed 90 cm/36"</i>), <i>tailored</i>	90 x 200 cm/36" x 80"	grey/white striped
1651	WendyLett, 80 cm/32" sliding surface (<i>bed 105 cm/42"</i>), <i>tailored</i>	105 x 200 cm/42" x 80"	grey/white striped

WendyLett2Way - drawsheet

(To be combined with WendyLett sliding sheet)

1620	WendyLett2Way	90 x 180 cm/36" x 72"	white
1631	WendyLett2Way	100 x 200 cm/40" x 80"	white/grey striped
1636	WendyLett2Way	140 x 200 cm/56" x 80"	white/grey striped

WendyLett4Way - sliding drawsheet

(To be combined with WendyLett sliding sheet)

1637	WendyLett4Way	140 x 200 cm/56" x 80"	grey/white checked
1647	WendyLett4Way	200 x 200 cm/80" x 80"	grey/white checked
1648	WendyLett4Way, inco	75 x 200 cm/29 x 80"	grey (inco) checked
1649	WendyLett4Way, inco with handles	140 x 120 cm/56" x 47"	grey (inco)
1657	WendyLett4Way, with handles	140 x 190 cm/56" x 75"	grey/white checked

WendyLean - pillowcase with handles

1680	WendyLean	50 x 60 cm/20" x 24"	grey/white striped
------	-----------	----------------------	--------------------

Contact your local distributor if you have any questions about the product and its use. See www.vancare.com for a complete list of distributors. Always make sure that you have the right version of the manual. The most recent editions of the manuals are available for downloading from our website, www.vancare.com

Manufactured for:

XPECTG: Tpe0

37373uUxgeVgn: 22/8; 6/6747Kvgtpgvly y y &ceptg&qo "

Cwtqt. P G'8: : 3: Hz<624/8; 6/5; ; 6G'o c n3phB xceptg&qo

WUC

VANCARE[®]
We move body and mind