


Vanderlift™ Components

- A. Model Name** Vanderlift™
- B. Four Point Floating Hanger Assembly** The 4 point self adjusting sling attachment hanger assembly is standard on Vancare lifts. A two point style is also available.
- C. Safety Latch** Prevents the loop of the sling from accidentally coming off.
- D. Actuator** Used to raise and lower the lift arm.
- E. Emergency Manual Lowering Device** Used to manually lower the lift in the event of power failure.
- F. Front Caster** Dual 3" Caster is hospital grade.
- G. Rear Caster** The heavy-duty 4" locking Caster is hospital grade.
- H. Auxiliary Up/Down Switch** Used to raise/lower the lift in the event of hand control failure.
- I. Base/Battery Box** Houses electronics for lift operation including the battery.
- J. Charging Cord** Hospital grade charging cord plugs into a hospital grade receptacle or into a standard wall outlet to charge the on-board batteries.
- K. Base Spread Handle** Used to open and close the legs of the lift.
- L. Leg Spread Position Bracket** Secures legs in three locking positions when spreading the lift legs.
- M. Handle Assembly** Handles used to maneuver the lift.
- N. Hand Control** Magnetized control used to raise and lower the lift.
- O. Emergency Stop Switch** Used to shut-off power to the lift.

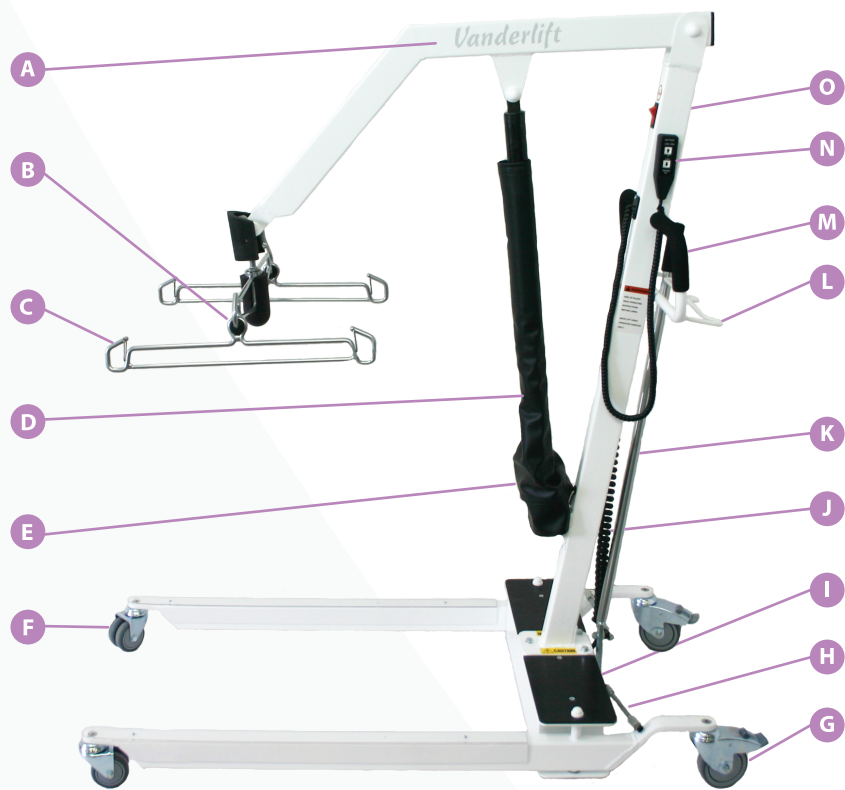


VANCARE®

We move body and mind

Vanderlift™

450, 600, 1000

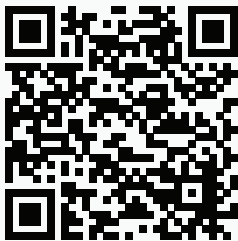

Scan the barcode above with a smart phone for more product information.

VANCARE, INC
1515 1st STREET
AURORA, NE 68818
TEL 800.694.4525
TEL 402.694.4525
FAX 402.694.3994
vancare.com



Made with Pride in the USA

You can trust Vancare

You can trust the Vanderlift™ to safely transfer those in your care, regardless of their range-of-motion limitations. Choose the model that best fits your needs:

- The Vanderlift 600 with a 600 lbs. capacity (page 1)
- The Vanderlift 1000 with a 1000 lbs. capacity (page 2)
- The Vanderlift 450 with a 450 lbs. capacity (back cover)

All models have a proven track record of reducing back injuries common among nursing staff - Minimizing lost work days, worker's compensation claims, staff turnover, and liability insurance premiums.



Check out the features that make the Vanderlift™ the brand of choice

- Magnetized hand-held control secures anywhere on the lift
- Auxiliary up/down switch provides alternative means of raising or lowering the resident/patient. (In case of a power failure, the resident/patient can be safely lowered with an emergency manual lowering system.)
- Onboard charger assures trouble free operation. The Vanderlift™ is simply plugged into a hospital grade or standard wall outlet to charge.
- Light blue, White Hammertone™, or Vancare Purple™ finish offer years of chip-resistant, non-peeling wear.
- Complies with all applicable federal and state regulations (i.e. OSHA, UL 2601-1, UL60601-1).
- Comes standard with self-adjusting floating four point sling attachment that offers ultimate comfort and security for the resident/patient.



The versatile Vanderlift™ Series Lift is the ultimate risk manager!

Whatever the transferring need - lifting a resident/patient from the floor or raising someone onto an elevated rehab table - the versatile Vanderlift™ Series Lift easily adapts. Even the largest and least ambulatory people can be transferred safely and comfortably.

World-class slings to meet your needs

Vancare's sling/belt selection is world class. With several hundred models, sizes and material choices, Vancare can provide the perfect fit for your patients and your budget. Vancare also proudly produces top quality "Made in the USA" products. When your patient's safety is at stake, you can trust a Vancare sling or belt.

