




## C-300 Specifications

- Lift Motor: 24 VDC
- Charger Input/Output: 100-240 VAC, 100 VAC:0.9 Amps / 24 VDC, 1.33 Amps
- Charger Model : SL POWER MENB1030A2403C01
- Batteries: NiMH —24 VDC (2 x 12 VDC) 1.8 AH-Custom Battery Pack.
- Lift Case: Flame Retardant ABS
- Hand Control: Digital( Operating force less than 1.12lbs or 5N)
- Lifting Range: Up to 1.96m
- Lift Weight: 5.7 lbs. (2.6 kg)
- Maximum Load: 300 lbs. (136 kg)
- Duty Cycle: 1 min "ON", 9 mins "OFF"
- Rated Performance: 20-21 lifts at 300 lbs.
- Max Sound Level: Raising Max load 56.1 dB , Lowering Max load 60.3 dB
- Service Life: 11000 Cycles or 5 years , whichever comes earlier.

## C-300 Dimensions

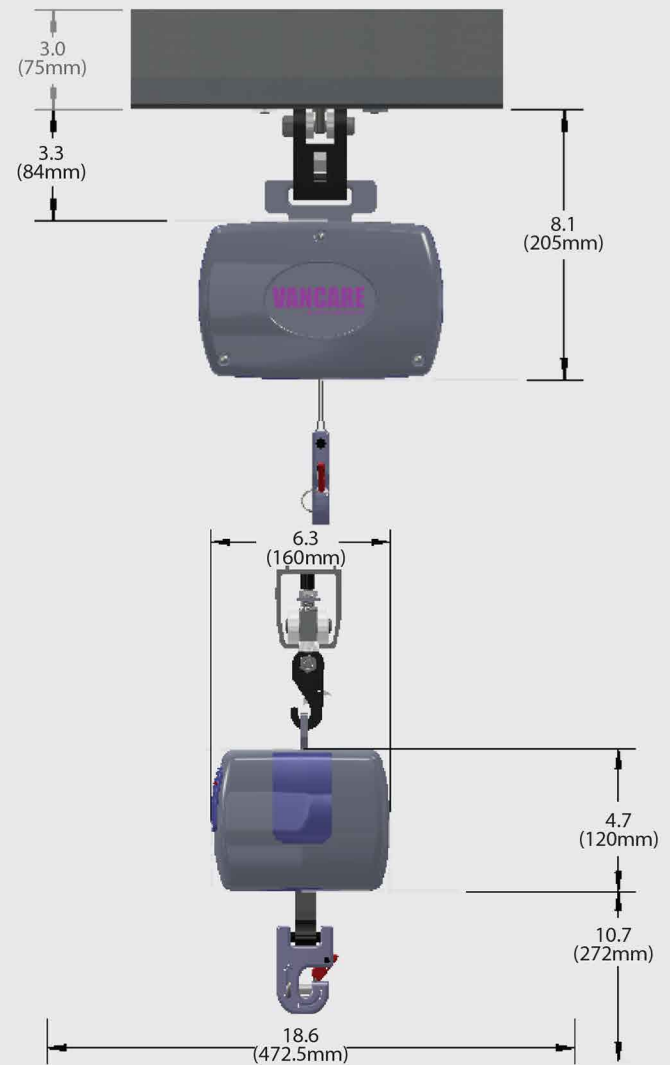

## C-300 Fixed/Movable Ceiling Lift 300 lb. Capacity

Committed to solving your safe patient handling challenges.

Vancare's newest extension of the C-Series design of ceiling lifts delivers the utmost in innovation, functionality and affordability. Ceiling lifts are widely regarded as the most effective transfer solution for both institutional and homecare settings. The C-300 ceiling lift has been designed and tested to the most rigorous industry standards. This product offers unparalleled value while compromising nothing with respect to its performance. With its 300lb. (136 kg) capacity and compact design, this lift will meet the needs of even the most demanding clinical environments. The C-300 is designed for use by a single caregiver offering quiet operation and effortless transfers. The Quick Release mechanism from the trolley furthers the simplicity and serviceability of this state of the art design.



The C-300 is equipped with a 3-button digital hand control with UP, DOWN and EMERGENCY DOWN buttons. The UP button is used to lift the patient and DOWN button is for lowering the patient.



Aurora, NE 68818  
Tel 800.694.4525  
Fax 402.694.3994  
Email [info@vancare.com](mailto:info@vancare.com)  
Web [vancare.com](http://vancare.com)

**VANCARE**<sup>®</sup>  
We move body and mind